





Cadenza CD-C18

30 x 4.6 mm
30 x 3 mm
30 x 2 mm


Application

No. TI021E


Agricultural Chemicals (1)

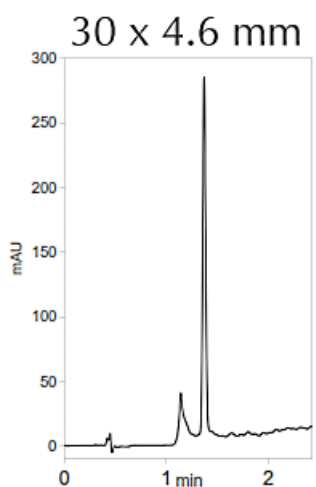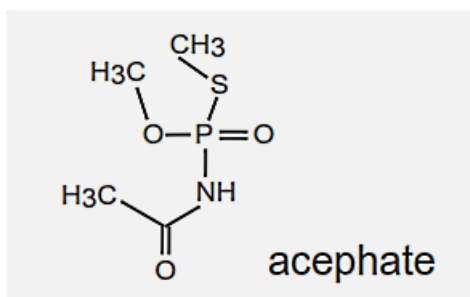

1.0 mL/min, 2.5M Pa
2.0 uL

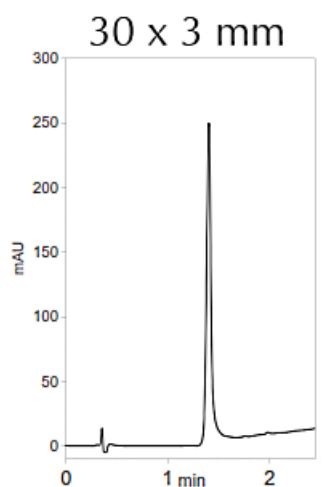

0.43 mL/min, 2.5 MPa
1.0 uL

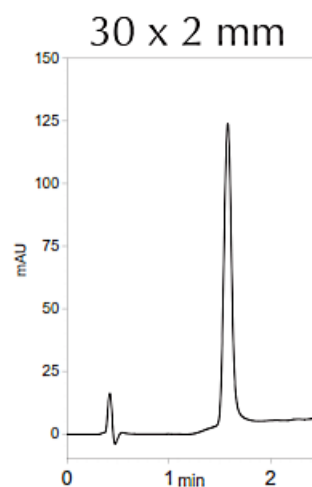

0.2 mL/min, 2.2 MPa
0.4 uL

Cadenza CD-C18
A: water / TFA = 100 / 0.1
B: acetonitrile / TFA = 100 / 0.1
0 - 20 %B (0-2 min)
UV at 220 nm, 37 °C