




Cadenza CL-C18
Unison UK-C18
Cadenza CD-C18
Unison UK-Silica

75 x 4.6 mm
250 x 4.6 mm


Application

No. TI198E


Antioxidants

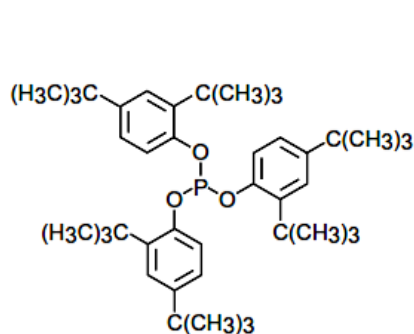


Irganox 1010

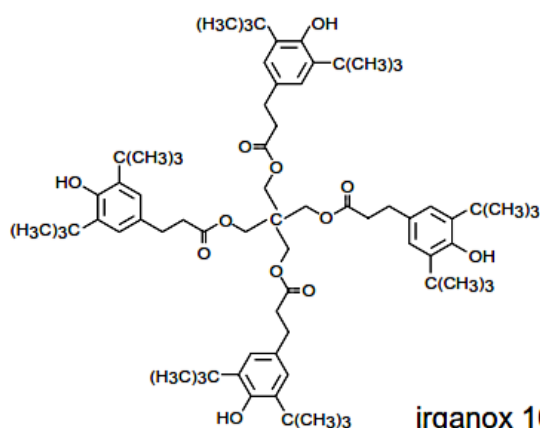

Irganox 168

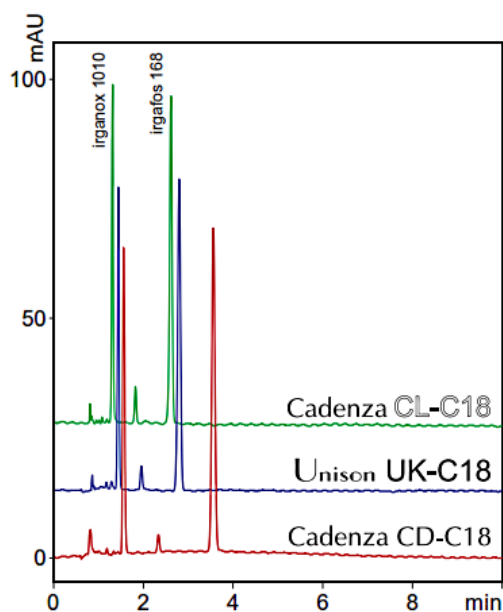

75 x 4.6 mm
acetonitrile / tetrahydrofuran
= 80 / 20
1 mL/min, 37 deg.C, 230nm

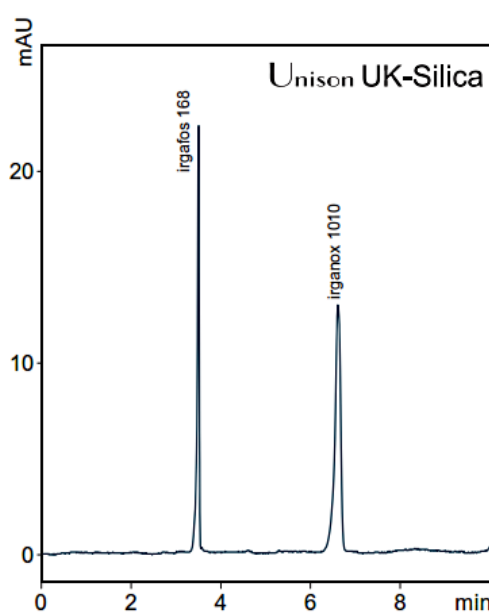

Unison UK-Silica
250 x 4.6 mm
hexane / ethyl acetate = 90 / 10
1 mL/min, 37 deg.C, 270nm